| | | |
|-------------------|---|----|
| •English | Addendum to the Instructions for Use | 1 |
| •Français | Addenda au manuel d'utilisation | 5 |
| •Deutsch | Nachtrag zur Gebrauchsanweisung | 9 |
| •Nederlands | Addendum op de gebruiksaanwijzing | 13 |
| •Español | Anexo a las Instrucciones de uso. | 17 |
| •Português | Adendo às Instruções de Uso. | 21 |
| •Italiano | Addendum alle Istruzioni d'uso | 25 |
| •Norsk | Tillegg til brukerhåndboken | 29 |
| •Svenska | Tillägg till bruksanvisningen. | 33 |
| •Suomi | Käyttöoppaan lisäys. | 37 |
| •Dansk | Tillæg til Brugerhåndbogen. | 41 |
| •Magyar | Kiegészítés a használati útmutatóhoz | 45 |
| •Polski | Dodatek do instrukcji obsługi | 49 |
| •Česky | Dodatek k Návodu k obsluze. | 53 |
| •Slovensky | Dodatek k Návodu k obsluze. | 57 |
| •Română | Anexă la instrucțiunile de utilizare | 61 |
| •Eestikeelne | Kasutusjuhendi lisa | 65 |
| •На русском языке | Приложение к инструкции по эксплуатации | 69 |
| •Қазақ тілінде | Пайдаланушы нұсқаулығына қосымша | 73 |
| •Ελληνικά | Συμπλήρωμα στις Οδηγίες Χρήσης | 77 |
| •Türkçe | Kullanım Talimatları'na Ek | 81 |
| •日本語 | ユーザーズ・ガイド補足. | 85 |
| •繁體中文 | 使用說明附註. | 89 |
| •简体中文 | 《使用说明》附录. | 93 |

Addendum

English

Addendum to the Instructions for Use

This addendum provides additional information for the Avalon FM20/FM30/FM40/FM50 Fetal Monitor Instructions for Use for F.01.xx (M2703-9001D). If your Instructions for Use have a different part number, you may discard this addendum, otherwise please store it with your monitor documentation.

OB TraceVue: via LAN

Both the monitor from which you are transferring a patient and the monitor to which you are transferring her must be switched on and connected to the OB TraceVue network during the patient transfer.

Physical Specifications

| Monitor Physical Specifications | M2704A/M2705A |
|---------------------------------|--|
| Water Ingress Protection Code | IP X1 (provided recorder paper drawer is shut) |

Monitor Performance Specifications

| Monitor Performance Specifications | M2702A/M2703A | M2704A/M2705A |
|------------------------------------|--|---------------|
| Global Speed | 6.25 mm/sec, 12.5 mm/sec, 25 mm/sec, 50 mm/sec | |

MECG Specifications

| ECG/Arrhythmia/ST Supplemental Information as required by IEC 60601 - 2- 27 | |
|---|--|
| Maximum T-wave amplitude | a = 8.0 mm |
| Heart Rate Averaging Method | Three different methods are used: Normally, heart rate is computed by averaging the 12 most recent RR intervals. For runs of PVCs, up to 8 RR intervals are averaged to compute the HR. If each of 3 consecutive RR intervals is greater than 1200 ms (that is, rate less than 50 bpm), then the 4 most recent RR intervals are averaged to compute the HR. |
| Display Update Rate | 2 s |
| Response Time of Heart Rate Meter to Change in Heart Rate | HR change from 80 to 120 bpm: Average: 9 seconds HR change from 80 to 40 bpm: Average: 10.5 seconds |
| Heart Rate Meter Accuracy and Response to Irregular Rhythm | Ventricular bigeminy: 63 bpm Slow alternating ventricular bigeminy: 46 bpm Rapid alternating ventricular bigeminy: 93 bpm Bidirectional systoles: 125 bpm |

Monitoring FHR Using DECG

CAUTION If a trace appears suspicious, verify the fetal heart rate with alternative means, e.g. auscultation, switch to fetal monitor's ultrasound transducer, sonography.

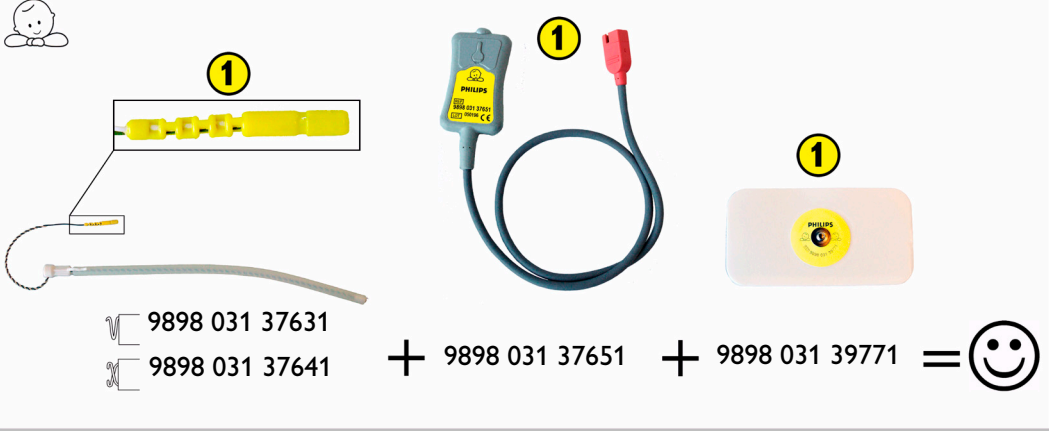
Troubleshooting

| Problem | Possible Causes | Solutions |
|--|--|--|
| DECG LEADS OFF displayed. Numeric is displayed with a - ? - ; INOP tone. See also Chapter 5, "Patient Alarms and INOPs". | Spiral electrode detached at connector. | Reconnect the fetal scalp electrode. |
| | Poor or no contact between leg attachment electrode and mother. | Check all connections. |
| | No contact between the DECG adapter cable and the leg attachment electrode. | Disconnect and reconnect the connector several times. |
| | No contact between the fetal scalp electrode connector and the DECG adapter cable. | Check all connections. Disconnect and reconnect the connector several times. If problem persists, use a new fetal scalp electrode. |
| Erratic trace. Erratic display. | No ECG signal. | Check for fetal demise. |
| | Poor contact between the reference electrode and the mother. | Use a new fetal scalp electrode if necessary. |
| Suspicious Pattern | Poor Contacts | Check all connections and electrodes. |
| | Interferences | Disconnect and reconnect the connector several times. |
| | Fetal arrhythmia. | Verify the fetal heart rate with alternative means |
| Signal quality indicator continuously shows a poor signal. | Fetal arrhythmia. | Ensure artifact suppression is off. |
| DECG EQUIP MALF displayed. | | See the chapter "Patient Alarms and INOPs". |
| DECG UNPLUGGED displayed. | | |
| DECG SIGNAL LOSS displayed. | | |

DECG Accessories: Component Compatibility

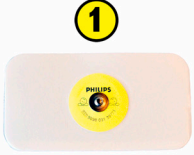
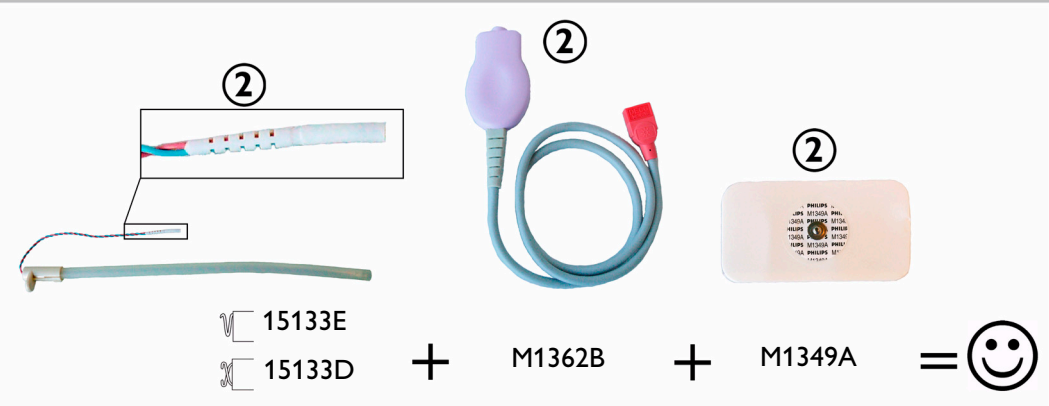
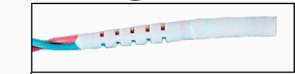
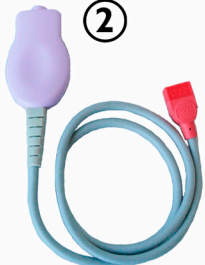
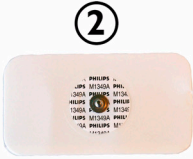
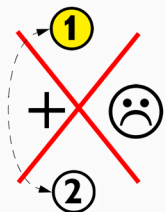
Use the following pictorial guide to check component compatibility for DECG accessories. **Do NOT** mix accessories from the New Philips DECG Solution (marked ①) with those from the QwikConnect Plus™ Solution (marked ②).

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